



Attention: Before operating BSB series clean booth, user must read this operating manual carefully, please! Keeping the manual for reference.

BSB Series Clean Booth

OPERATING MANUAL

ISO9001: 2000 CERTIFIED

SUZHOU CLEANING TECH. RESEARCH INSTITUTE

BAISHEN TECHNOLOGY (SUZHOU) CO., LTD.

Welcome to use BSB series clean booth, thank you warmly entering into our company's consumer team!

1.General

Along with science and technology's developing rapidly, air clean technology has already been used extensively in production and scientific research department of electron, spaceflight, aviation, precision instrument, apparatus & meter, chemical industry, light industry, medicine and foodstuff.

The Clean booth is a kind of vertical unidirectional airflow cleaning device, It offers great effects on improving technical conditions and raising the accuracy purity and reality of the products. Along with science technical rapidly develop, air clean technology already cover extensively should be used for laboratory, electron technology, aviation and space, precision instrument.

2 Technical specifications

1) Cleanliness :100 class (@ $\geq 0.5 \mu m$ FS 209E)

2) Mean falling bacteria : ≤ 0.5 piece / culture container. hour (ϕ 90mm culture container)

3) Mean flow rate (m/s): 0.3~0.6 (adjustable)

4) Noise level ≤ 65 dB (A)

5) Vibration/half peak (μm) : $\leq 5 \mu m$

6) Illumination (Lx) ≥ 300 Lx

7) Power supply: 220V/50Hz

Model	BSB-880	BSB-1360	BSB-1980
Outline Dimension (mm)	1040 ×700×1500	1480 ×700×1500	2140 ×700×2150
Internal Dimension (mm)	880 ×700×1150	1320 ×700×1500	1980 ×700×1800
Lluorometry	40W×①	30W×②	30W×②
Size & No. of HEPA (mm)	820 ×610×①	610×610×②	610×610×③
Weight	80Kg	100Kg	180Kg
MAX. consumption (VA)	400	600	900

Note: Thick of filter can be selected 50 or 69mm.

3 Operating Principle

BSB series belongs to vertical unidirectional airflow, Possess positive and negative press wind channels circular system:

under the usual condition(positive press) After being firstly filtered by the roughing air filter the indoor air is sucked into the static pressure case through variable velocity centrifugal machine. It is secondly filtered by the high efficiency particulate air filter (HEPA) to reach the bench work surface, The cleaning flow through the outlet of the HEPA with a mean and certain crosscut velocity through the operating range takes the dust particles and biological particles away to make an operating range with no bacteria.

4.Installation

- 4.1 His device must be installed in the cleaning room (100000class at lowest)
- 4.2 Please confirm that this product's rated voltage, frequency conform to power supply's parameter, prevent electric element and device from damage and lead to fire.
- 4.3 The power cord must be connecting to the socket which is grounded.

5.Usage

- 5.1 Turn-on the power supply, open the ultraviolet bacteria-killer lamp, after 30 minutes, close ultraviolet bacteria-killer lamp and open the fluorescent lamp. Then, switch on the flow motor and start operating after 20 minutes, open the fluorescent lamp when operating
- 5.2 If the working cabin has not been used for a long time, before operating its circle must be cleaned with super vacuum cleaners or other tools, which cannot make fiber. At last bactericidal agent or ultraviolet rays should be used to kill the bacteria.
- 5.3 Unnecessary things are not allowed to be put on the working table to assure the cleaning flow track in the working range not to be disturbed.
- 5.4 In the cleaning operating range any obvious action to disturb the flow track is not

allowed.

5.5 Operating area's usage temperature can not be more than 60 ℃ .

6.Maintenance


6.1 According to environmental clean extent, the coarse filter cloth should be regularly (generally every 2-3months) cleaned by removing (Dacron non-woven cloth) or replacing.

6.2 it must be replaced HEPA when the high-speed range of fan could not make working area wind speed reach to 0.30m/s,

6.3 when replacing fuse, ensure capacity same as original to avoid fire.

4.when checking and maintenance, please cut off power switch and pull out the power plug.

7 Operation panel function illustrate

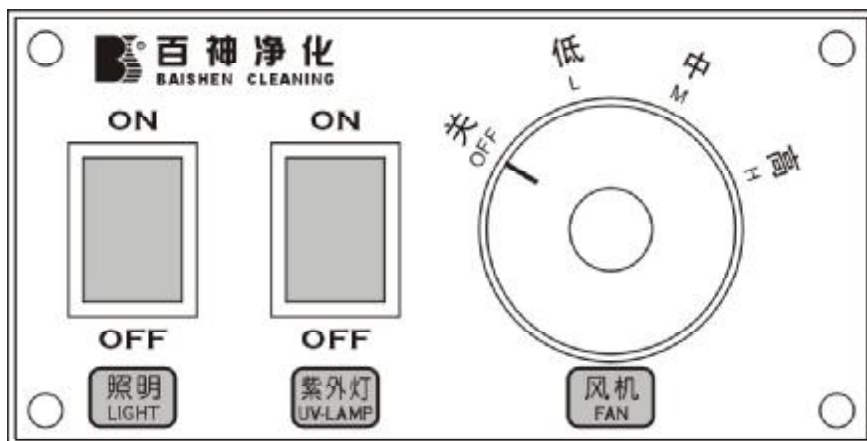

Remark: if it has dimension deflection, according to upgrade version.

BSB Series Clean Bench OPERATING MANUAL

**SUZHOU CLEANING TECH. RESEARCH INSTITUTE
BAISHEN TECHNOLOGY (SUZHOU) CO., LTD.**

Add: 31 Qihong Road, Tongyuan Road ,SIP Suzhou

Zip: 215006

Phone: +86-512-62522538

Fax: +86-512-62522938

[http: //www.baishen.com](http://www.baishen.com)

E-mail: clean@baishen.com

2004 Printed in Suzhou